PDA Training Course Extractables & Leachables 25-26 April 2024

The AET Challenge for Large Volume Parenterials (LVPs): Extractables Simulation Studies and How to Design Them

Dr. Dennis Jenke









Training Outline

- What is an Large Volume Parenteral (LVP)?
- The LVP Challenge How Low Can You Go?
- The Simulation Study as a Means of Addressing the LVP Challenge
- Design Parameters for Effective Simulation Studies
 - Extraction Solvent Composition
 - Temperature and Duration
 - Stoichiometry





Pharmaceutical Dosage Forms

Different Dosage Forms








Parenteral Dosage Forms

Parenteral drug products are injected through the skin or other external boundary tissue, or implanted within the body, to allow the direct administration of the active drug substance(s) into blood vessels, organs, tissues, or lesions. Parenteral dosage forms include solutions, suspensions, emulsions, sterile powders for solutions and suspensions (including liposomes), implants (including microparticles), and products that consist of both a drug and a device such as drug-eluting stents.

Routes of Parenteral Administration

- Intravenous injections and infusions
- Subcutaneous injections
- Intramuscular injections
- Intradermal injections
- Intra-arterial injections
- Intra-cardic injections
- Intraspinal injections
- Intra-articular injections

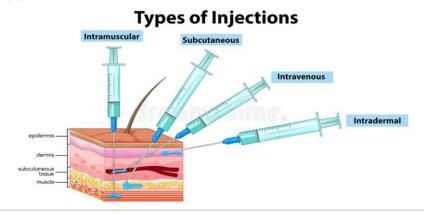





What is an LVP?

A **single-dose injection** that is intended **for intravenous use** and is packaged in containers labeled as containing **more than 100 mL**.


Characteristics of LVPs

- Packaged in glass bottles or in large volume flexible containers.
- May contain greater than 100 mL to greater than 1 or 2 L
- Sterile (e.g., many LVP are sterilized in their container via heat, although some are sterile-filled)
- Pyrogen-Free
- Essentially free of particulate matter
- No anti-microbial agents
- Isotonicity
- Longer term use





Types and Compositions of LVPs

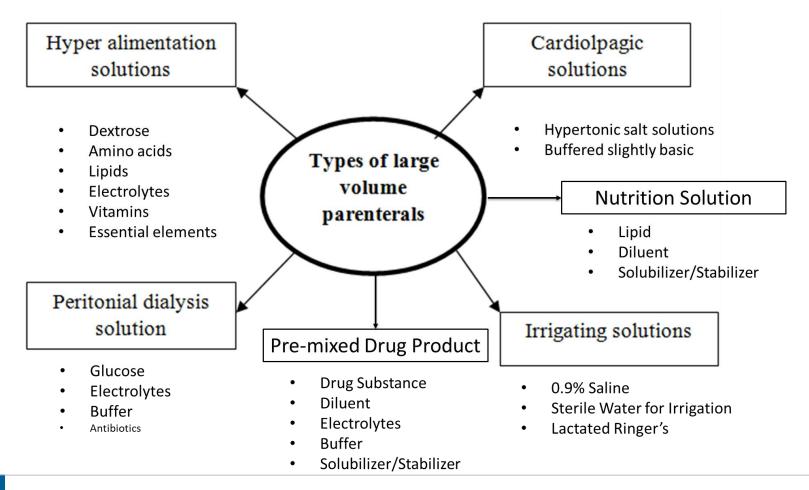



Generic Composition of an LVP

Vehicle

- Water
- Water miscible vehicle
- Non-aqueous vehicle
- Solid Vehicle

Active ingredient

Added substances

- Tonicity Adjusters
 - Electrolyte, NaCl, 0.5 0.9%
 - Non-electrolyte, Dextrose, 4 5%
- Buffers
 - Acetate/Citrate, pH 3 6
 - Phosphate, pH 6 8
 - Glutamate, pH 8 10
- Antioxidant, 0.1 0.5 %
- Preservatives, 1 2%
- Complexing agents, 0.01 0.05%
- Surfactants, 0.05 0.5%
- Competitive Binders, variable
- Antimicrobial agents, 0.01%
- Cryoprotectors/Lyoprotectors (Bulking agents), e.g., 1 10%
- Etc.





Packaging for LVPs



Container

Printing (label)

Access Port (elastomeric septum not shown)





Auxiliary
Items
(e.g. oxygen absorber)





The LVP Challenge for Leachables

Among the numerous characteristics that differentiate Large Volume Parenterals (LVPs) from other dosage forms, their **composition** and **large daily dose volume** are particularly noteworthy because of the practical implications of composition and dose volume to the safety assessment of packaging system leachables.









versus





The LVP Challenge for Leachables

Composition










Daily Dose Volume

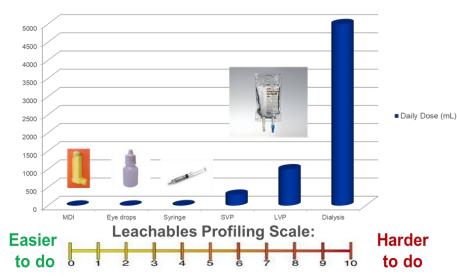






Why is Daily Dose Volume Important?

One of the most basic concepts in toxicological assessment is that:

"The dose makes the poison."

A substance can adversely affect health only if the amount of the substance to which an individual is exposed (dose) exceeds a tolerable threshold.

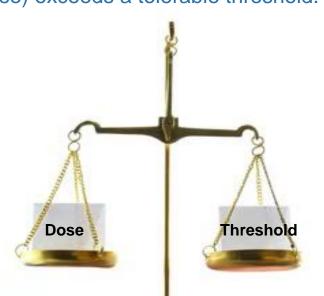

Paracelsus, the "Father" of modern toxicology


Dose = concentration in medication x volume of medication used





Why is Daily Dose Volume Important?

The value of the Safety Threshold decreases in direct proportion to the increase in Daily Dose Volume.

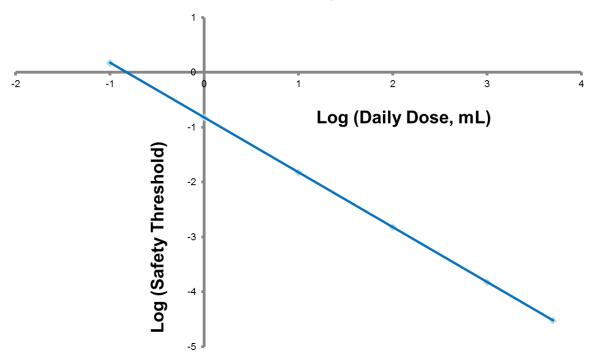






Daily Dose Volume and the AET



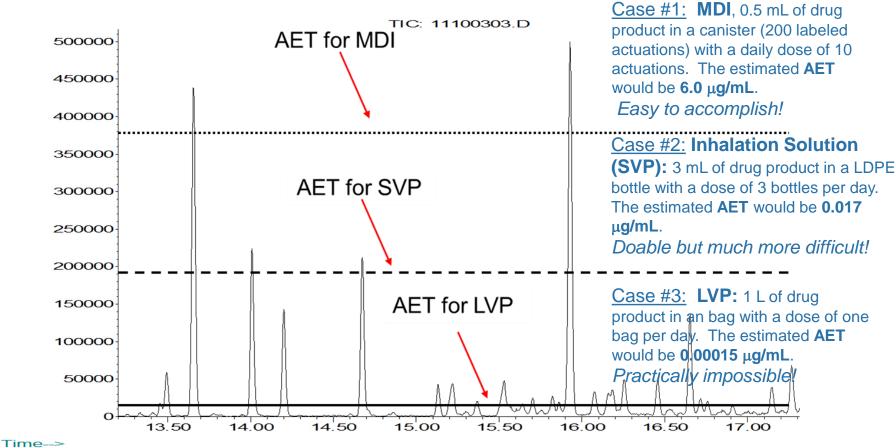



pda.org



The LVP Challenge: How Low Can You Go?



AETs for LVPs may be so low that even state of the art, best demonstrated practice analytical methods may not be able to accomplish the functions of discovery and identification for all necessary leachables.

If leachables cannot be detected and identified then obviously they cannot be toxicologically assessed by numerical means and thus their potential safety impact cannot be established by such numerical means.





The Simulation Study Goes Low!

1. The drug product formulation has been replaced with one or more simulating solvents that are easier to test.







2. The actual use conditions of contact have been accelerated.



3. The test article may have been altered (somewhat) to provide an exaggerated and presumably worst case.












Simulation Study Objectives

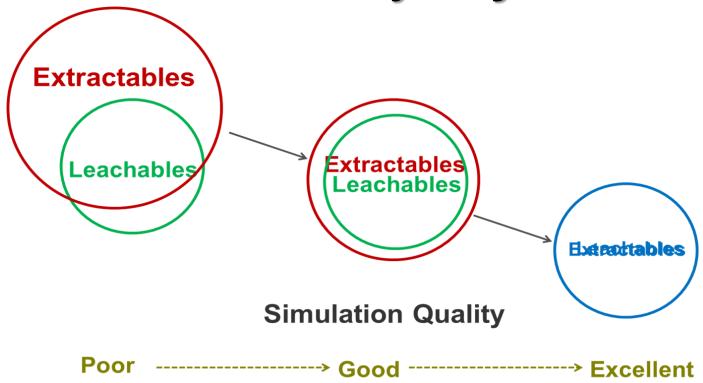


An extractables profile that is the same as a drug product's leachables profile (meaning that the extractables profile includes all the members of the leachables profile with extractables levels being greater than or equal to the leachables levels).





Key Design Parameters to Simulate

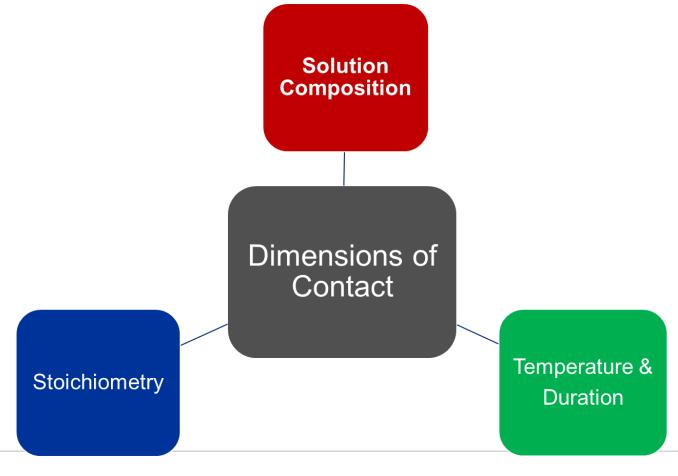





Simulating Solution Composition

Solution Composition

1. Polarity

2. pH

3. "Reactivity"





Simulating Polarity

- A leachable will accumulate in a drug product to a level dictated by its solubility in the drug product.
- A leachable's solubility in a drug product will depend on the "polarity" of the leachable and the drug product ("Like dissolves like").


Means of Establishing a Solution's Polarity:

- 1. Polarity Tables for Solvents
- 2. Correlation with Measurable Fundamental Properties Dielectric Constant
- 3. Use of Polarity Markers (e.g., solvatochromic Reichardt's dye)
- 4. Experimental Determination via "Extraction Power" Scales





An "Extraction Power" Scale

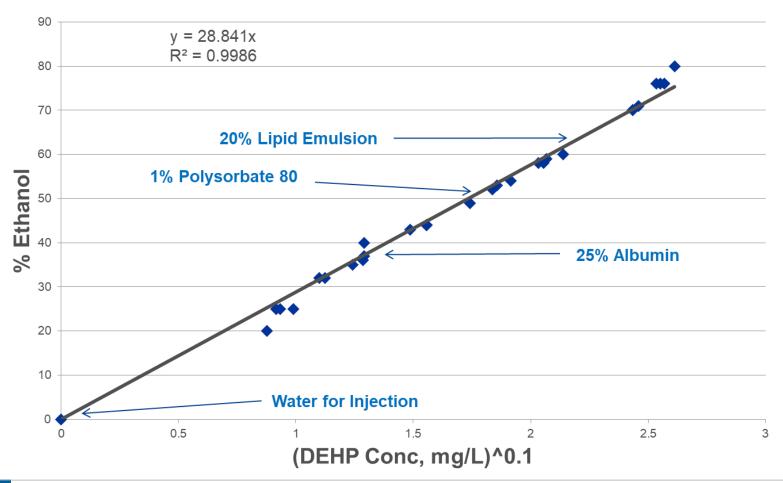





Realistic Expectations About Simulating Polarity

Realistically,


- Many LVP's are aqueous solutions containing inorganic and organic salts, acids, bases and organic compounds of high aqueous solubility (e.g., dextrose) and no solubility enhancers (e.g., Tween 80). In these cases, use of polar simulating extraction solvents is appropriate.
- A few LVPs, for example nutritional products, may contain lipids and proteins in relatively high concentrations (e.g., up to 20 - 25% by weight). As the presence of these lipids and proteins makes the products more non-polar, use of semi-polar extraction solvents (such as alcohol/water mixtures) may be appropriate.
- LVP's are NEVER non-polar solvents and thus non-polar extraction solvents are NEVER appropriate.





Addressing Solution pH

- A leachable will accumulate in a drug product to a level dictated by its solubility in the drug product.
- The solubility of an acidic or basic leachable will depend on the acid/base dissociation constant (pK_a) of the leachable and the pH of the drug product.

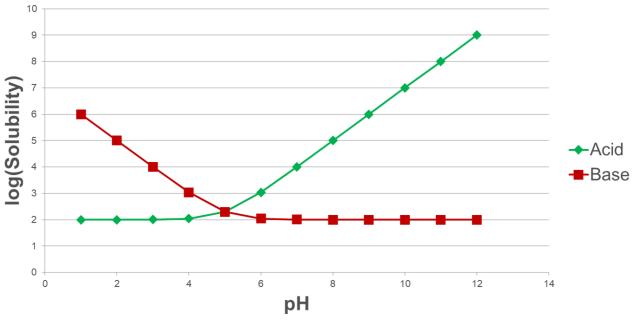





Simulating pH

- Most common acidic leachables have a pK_a of 7 or less.
- Most common basic leachables have a pK_a of 3 or more.
- Most aqueous drug products have a pH between 3 and 9.

Therefore:

- <u>Two simulating solvents</u>, one prepared at <u>pH 3</u> and one prepared at <u>pH 10</u>, reasonably bracket the universe of leachables and drug products, although acceptions may require more extreme pH values.
- There is little or no value in using a simulating solvent with an intermediate pH.
- If your drug product(s) have a more narrow pH range, then use a range that is more appropriate.





Accelerating Shelf-life: Temperature and Duration

Temperature and Duration






Accelerating an Extraction

The higher the temperature, the longer the extraction and the larger the diffusion coefficient:

- 1. The larger will be the leachable's concentration in the drug product.
- 2. The more likely an equilibrium leachable concentration will be achieved.

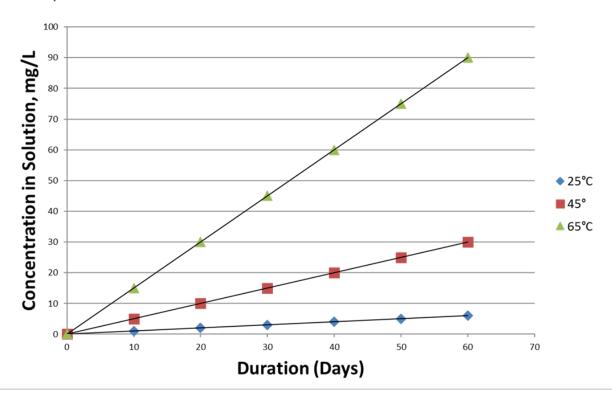





Estimating Accelerated Conditions (1)

 ASTM F1980-16: Standard Guide for Accelerated Aging of Sterile Barrier Systems for Medical Devices.

Accelerated Aging Time at T2 = Actual Aging Time at T1 ÷ C

$$C = Q_{10}^{[(T2-T1)/10]}$$

where $Q_{10} = 10^{\circ}C$ Reaction Rate Constant

T2 = accelerating temperature (°C)

T1 = actual temperature of contact (°C)





Estimating Accelerated Conditions (2)

2. "Factor 10 Rule"¹: This factor 10 rule is based on the observation that activation energies for migrating substances in polymers relevant to packaging are typically in the range of 80 to 100 kJ/mole. In such a circumstance, the diffusion coefficient increases by roughly an order of magnitude for every 20°C increase in contact temperature. Thus for example, the migration rate at 40°C is ten times faster than the migration rate at 20°C.

Accelerated Aging Time at T2 = Actual Aging Time at T1 ÷ C


 $C = 10^{[(T2 - T1)/20]}$





Which Approach to Use?

Acceleration of a Two-Year (730 days) Ambient Temperature Shelf-life

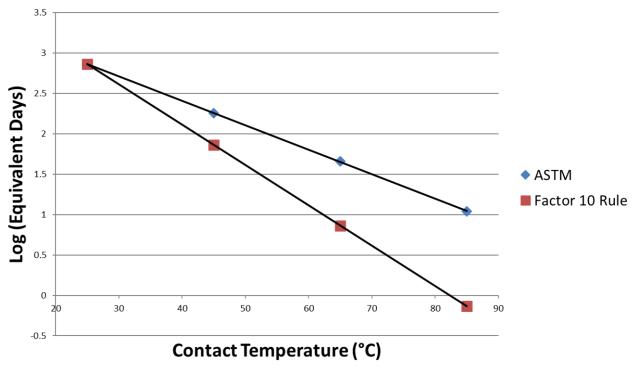

The ASTM approach produces the longest duration is thus is the most conservative.





Example Acceleration Calculation

The time (t_2) required for an extractable to reach a certain concentration at a temperature T2 can be estimated from the time (t_1) required for the same extractable to reach the same concentration at a reference temperature T1 using the following equation, although exceptions will occur:

$$\mathbf{t}_2 = \mathbf{t}_1 \div \mathbf{10}^{[(T2-T1)/20]}$$

For example, if the time it takes for an extractable to achieve a concentration of 2.0 mg/L at 25°C is 10 hours, the time it takes for the same extractable to achieved the same concentration of 2.0 mg/L at 45°C will be:

$$t_{45} = t_{25} \div 10^{[(45-25)/20]}$$

 $t_{45} = 10 \text{ hours} \div 10^1$



$$t_{45} = 1$$
 hour



Extraction Stoichiometry

Stoichiometry

- 1. Surface area/Solution volume
- 2. Material weight/Solution volume





Stoichiometry Fallacies:

Its all about surface area.

In fact, what is generally attributed to surface area effects is actually due to changes in the amount (mass) of the extracted item.

2. As the surface area to solution volume ratio increases, the concentration of extractables will increase in the same linear and 1 to 1 manner for all extractables.

In fact, the relationship between the ratio and the concentration depends on the plastic/solution partition coefficient of the extractables in question.

- "Solution-loving" extractables' concentrations will increase in proportion to increased ratio (but not 1 to 1)
- "Plastic-loving" extractables' concentration will change very little as the ratio increases.





The Stoichiometry Equation

$$C_{l,e} = m_{l,e}/V_l = m_{p,o}/[V_l + (k_{p/l} \times SA_p \times t_p)]$$

Where C is the extractable's concentration,


- m is the mass of the extractable in either phase,
- SA is the surface area of the sample being extracted,
- · t is the thickness of the sample being extracted,
- k_{p/l} is the extractable's plastic/solution partition coefficient,
- V is the volume of either phase, and
- the subscripts p, l, e and o refer to the plastic phase, the liquid phase, equilibrium and original respectively





The Real Effect of SA/V:

Theoretical Relationship between the Material Surface Area and the Concentration of an Extractable in an Extracting Solution at a Constant Extracting Solution Volume.

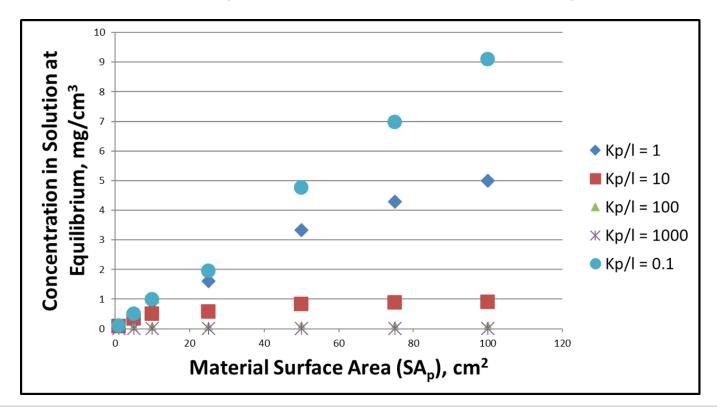





This is <u>not</u> a Simulated Extraction!



Once the scaling factor of a simulated extraction versus commercial drug product exceeds 3, then the effect of the partition coefficient is so great that the scaling factor cannot be used to extrapolate simulated extraction data to the commercial drug product.





In Review

- An LVP is an injected dosage form with a high unit volume and frequently complicated composition.
- A properly designed and implemented extractables simulation study produces an extractables profile
 that is equal to or slightly exaggerated than the leachables profile for a packaged drug product.
- Critical design parameters for a simulation study include:
 - Solution Composition
 - Temperature and Duration
 - Stoichiometry
- In considering Solution Composition, the aspects of "polarity", pH and "reactivity" should be considered. Of these three, "polarity" and pH are relatively straightforward, while "reactivity" needs further consideration.
- In considering Temperature and Duration, certain mathematical conventions can be quite useful in terms of accelerating leaching.
- In considering Stoichiometry, it is noted that in many cases the surface area to solution volume ratio is just another way of saying material weight to solution volume. More importantly, the assumption of a linear relationship between stoichiometry and leachables accumulation may or may not be true.





Key References

- <1664> Assessment of Drug Product Leachables Associated with Pharmaceutical Packaging/Delivery Systems. USP 38 – NF 33 (First Supplement), pp. 7181 – 7193. August 1, 2015.
- 2. Jenke, D.; Liu, N.; Hua, Y.; Swanson, S.; Bogseth, R. A means of establishing and justifying binary ethanol/water mixtures as simulating solvents in extractables studies. *PDA J Pharm Sci Technol.* 69(3): 366-38 2(2015).
- 3. Jenke, D. Establishing the proper pH of simulating solvents used in organic extractables assessments for packaging systems and their materials of construction used with aqueous parenteral drug products. *Pharm Outsourcing.* 15(4):20, 22, 24-27 (2014).
- 4. ASTM F1980-07 (Reapproved 2016): Standard Guide for Accelerated Aging of Sterile Barrier Systems for Medical Devices.
- 5. R. Franz, A. Stormer. Migration of Plastic Constituents. In Plastic Packaging: Interactions with Foods and Pharmaceuticals. Wiley-VCH; Second Edition, 2008, pp. 368.
- 6. R. Franz, A. Stormer. Migration of Plastic Constituents. In Plastic Packaging: Interactions with Foods and Pharmaceuticals. Wiley-VCH; Second Edition, 2008, pp. 370.
- 7. Jenke, D; Rabinow, B. Proper accounting for surface area to solution volume ratios in exaggerated extractions. *PDA J Pharm Sci Technol.* 71(3): 225-233 (2017).
- 8. Jenke, D. Application of Arrhenius Kinetics to Acceleration of Controlled Extraction Studies. *PDA J Pharm Sci Technol.* 73(2): 135-168 (2019).





Q&A

Thank you!



Contact the presenter at:

dennisjenke@triadscientificsolutions.com www.triadscientificsolutions.com



